

# HOTEL-STYLE RECLINERS FOR REST & CARE



Our range of REST AND HEALTHCARE RECLINERS is adapted for the Day Surgery, Day Care, General Medicine, Oncology and Dialysis units. Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.



## CONTENTS

▶ CAREXIA FP

▶ CAREXIA FPE

▶ CAREXIA FPU

▶ CAREXIA FPU E

▶ TIM

▶ VELVET

▶ SWANN

▶ COLOR CHART

Comfortable

Easy to use

Adaptable Various options available

Versatile Adapted to different units



## CAREXIA FP FIXED HEIGHT TREATMENT CHAIR

Carexia FP, is a fixed height versatile treatment chair.

Its soft and welcoming shapes bring comfort, relaxation and serenity to the patient.

Its pure and harmonious design is the vector of a positive image for your unit.




# CAREXIA FP

## CAREXIA FP DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	76.4 cm (folded) 170 cm (open)	69.5 cm 90 cm (with FPOAC3)	130 cm
SEAT	50 cm	57 cm	55 cm on seat
INCLINE FOLDER	85 to 60°		

SAFE WORKING LOAD: 200 kg

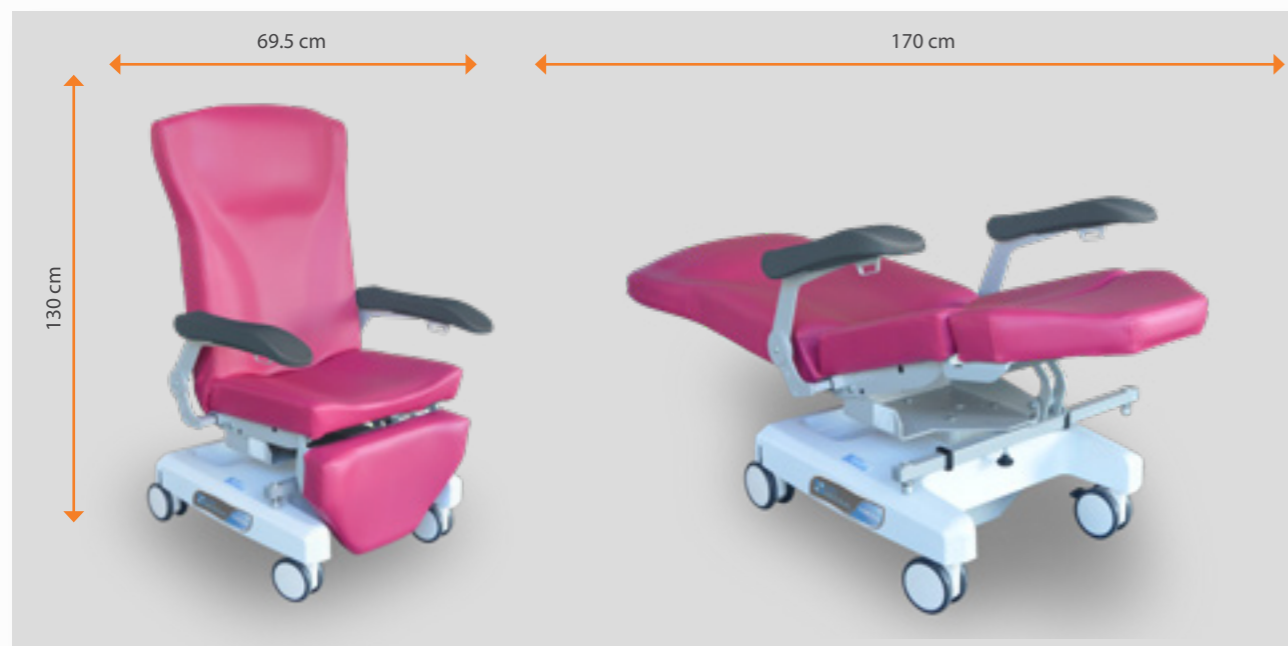

NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, General Medicine, Oncology and Dialysis units.


Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Versatile treatment chair.

- Heavy duty design. Steel frame with powder coating. ABS cover.
- Fixed height : 55 cm.
- A Ergonomic shaped backrest and seat. PVC coated jersey.
- High resilience foam 35 kg/m3.
- Integrated headrest.
- Backrest adjustable from 85 to 60° with gas spring assistance (85 to 6° with legrest option)
- Trend position with optional legrest.
- 4 self lockable casters with double running wheels
- ABS basket
- Polyurethan armrests, can be folded up




4 self lockable casters with double running wheels

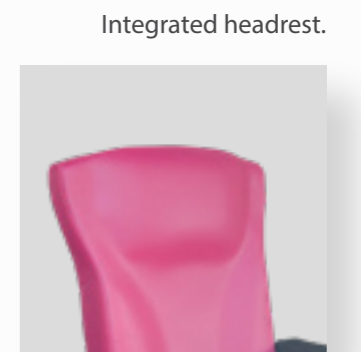

Integrated headrest.



ABS basket

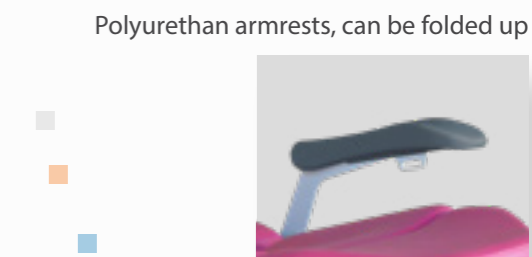


Polyurethan armrests, can be folded up

Faux leather covers : URBAN TEXTILE



ABSINTHE TOMATO

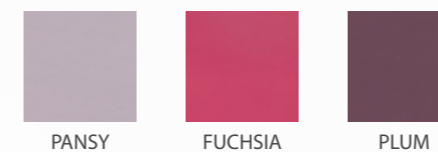

PANSY FUCHSIA PLUM




DEEP OCEAN PETROL GREY

URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %  
Thickness: 1,2 mm +/- 0,1 mm  
Effective against bacteria (MRSA)  
Phthalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof



# ACCESSORIES

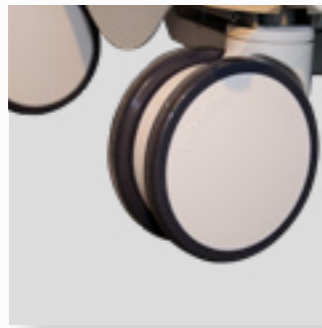

4 castors Ø 125 mm with central lock  
FPORS1



4 double runing wheel castors Ø 150 mm with central lock  
FPORS2




Push bar & Paper roll holder  
FPOBP & FPOSR



Short legrest synchronized with the backrest. Trend position.  
FPORJ



2 holders for IV pole (necessary for TSR2)  
FPOSA

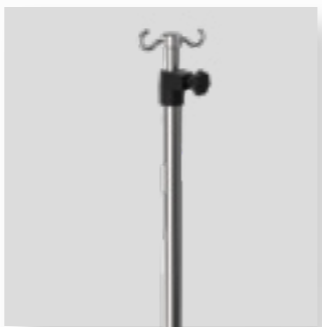

Removable IV pole, 2-stage, 2 hooks  
TSR2



Comfort pillow (grey)  
ACC-CTG



Armrest comfort (pair)  
FPOAC4



Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)  
FPOAC3

REFERENCE	DESCRIPTION
FPORS1	4 castors Ø 125 mm with central lock
FPORS2	4 double runing wheel castors Ø 150 mm with central lock
FPOBP	Push bar
FPORJ	Short legrest synchronized with the backrest. Trend position.
FPOSA	2 holders for IV pole (necessary for TSR2)
TSR2	Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)
TSR4	Removable IV pole, 2-stage, 4 hooks (FPOSA mandatory)
ACC-CTG	Comfort pillow (grey)
FPOAC4	Armrest comfort (pair)
FPOAC3	Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)
FPOSR	Paper roll holder
FPOGC	High comfort upholstery with 2 cm thickness memory foam
SUPTAB	Fixed I-Pad support with antitheft device

**Comfortable**

**Easy to use** Electric movements

**Adaptable** Various options available

**Versatile** Adapted to different units



## CAREXIA FPE FIXED HEIGHT TREATMENT CHAIR

CAREXIA is a versatile treatment chair.

Its soft and welcoming shapes bring comfort, relaxation and serenity to the patient.

Pure and harmonious design is the vector of a positive image for your unit.

CAREXIA is a chair designed for post-surgery rest, blood sampling, care or examinations, chemotherapy or hemodialysis.




# CAREXIA FPE

## CAREXIA FPE DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	76.4 cm (folded) 170 cm (open)	69.5 cm 90 cm (with FPOAC3)	130 cm
SEAT	50 cm	57 cm	55 cm on seat
INCLINE FOLDER	85 to 6°		

SAFE WORKING LOAD: 200 kg

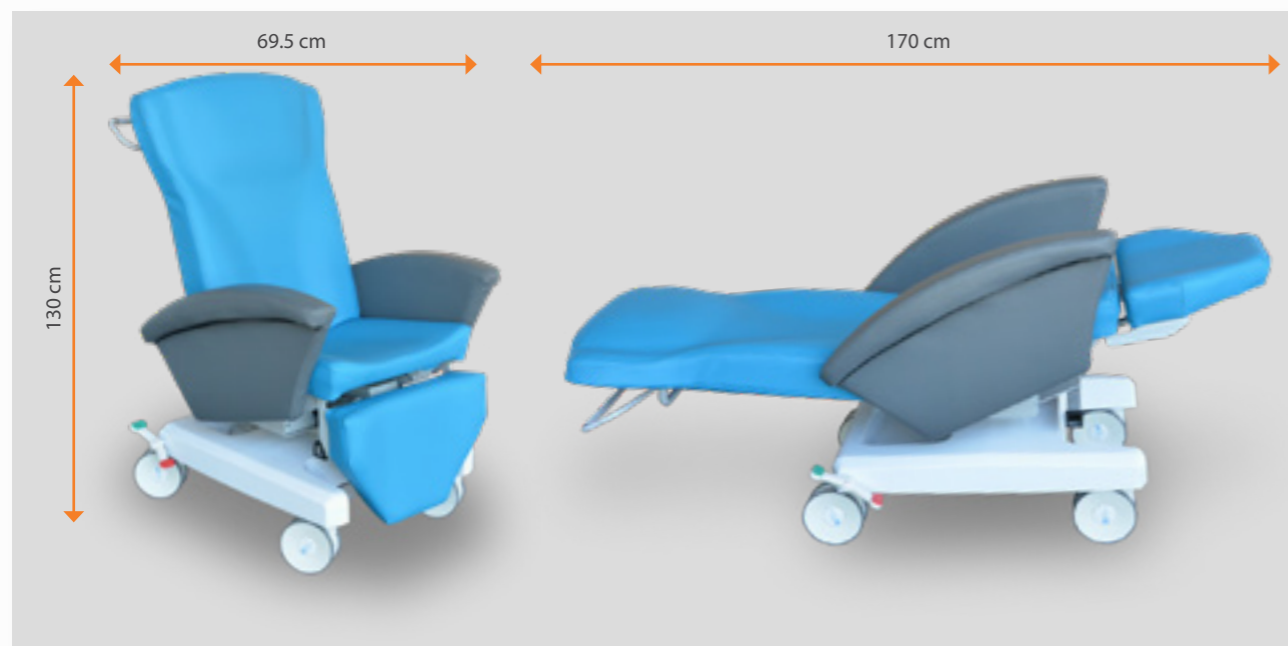


NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Multi-purpose treatment chair with electric movements.

- Heavy duty design. Steel frame with powder coating. ABS cover
- Fixed height : 55 cm
- Ergonomic shaped backrest and seat. PVC coated jersey
- High resilience foam 35 kg/m<sup>3</sup>
- Integrated headrest
- Electric backrest and legrest
- Trend position
- 4 self lockable casters with double running wheels
- ABS basket
- Polyurethan armrests, can be folded up

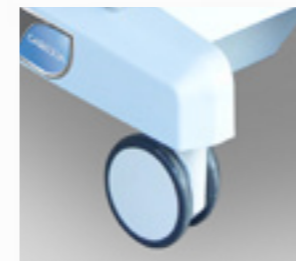


4 self lockable casters with double running wheels



ABS basket



Integrated headrest

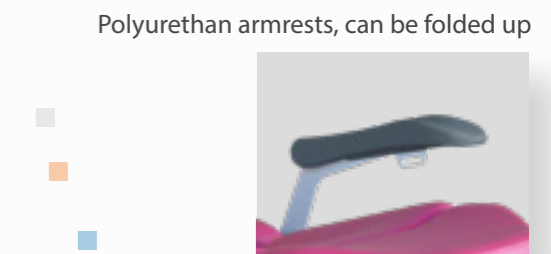


Polyurethan armrests, can be folded up

Faux leather covers : URBAN TEXTILE



ABSINTHE TOMATO

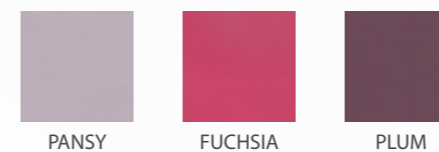

PANSY FUCHSIA PLUM



DEEP OCEAN PETROL GREY

URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %  
Thickness: 1,2 mm +/- 0,1 mm  
Effective against bacteria (MRSA)  
Phthalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof



# CAREXIA FPE




## ACCESSORIES

REFERENCE	DESCRIPTION
FPORS1	4 castors Ø 125 mm with central lock
FPORS2	4 double runing wheel castors Ø 150 mm with central lock
FPOBP	Push bar
FPORJL	Long legrest
FPORP	Electric adjustable footrest (only for long legrest) (incompatible with FPOAC4)
FPOCP1	Foot pedal control for relax and trend position
FPOSA	2 holders for IV pole (necessary for TSR2)
TSR2	Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)
TSR4	Removable IV pole, 2-stage, 4 hooks (FPOSA mandatory)
ACC-CTG	Comfort pillow (grey)
FPOAC4	Armrest comfort (pair)
FPOAC3	Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)
FPOSR	Paper roll holder
FPOGC	High comfort upholstery with 2 cm thickness memory foam
FPOLAMP	Reading lamp DEL
FPOBA	Battery 2 Ah
SUP-USB	USB Port
PROTEC RJL	Removable Crystal protection for footrest
PROTEC RP	Removable Crystal protection for footrest
SUPTAB	Fixed I-Pad support with antitheft device



# ACCESSORIES

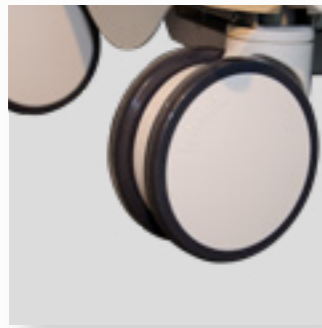

4 castors Ø 125 mm with central lock  
FPORS1



4 double running wheel castors Ø 150 mm with central lock  
FPORS2



Push bar & Paper roll holder  
FPOBP & FPOSR




Long legrest & Electric adjustable footrest (only for long legrest)  
FPORJL & FPORP



Foot pedal control for relax and trend position  
FPOCP1



2 holders for IV pole (necessary for TSR2)  
FPOSA

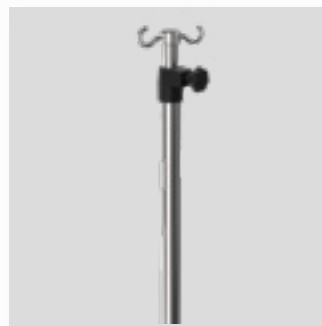

Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)  
TSR2



Comfort pillow (grey)  
ACC-CTG



Armrest comfort (pair)  
FPOAC4




Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair) FPOAC3



Reading lamp DEL  
FPOLAMP



USB Port  
SUP-USB

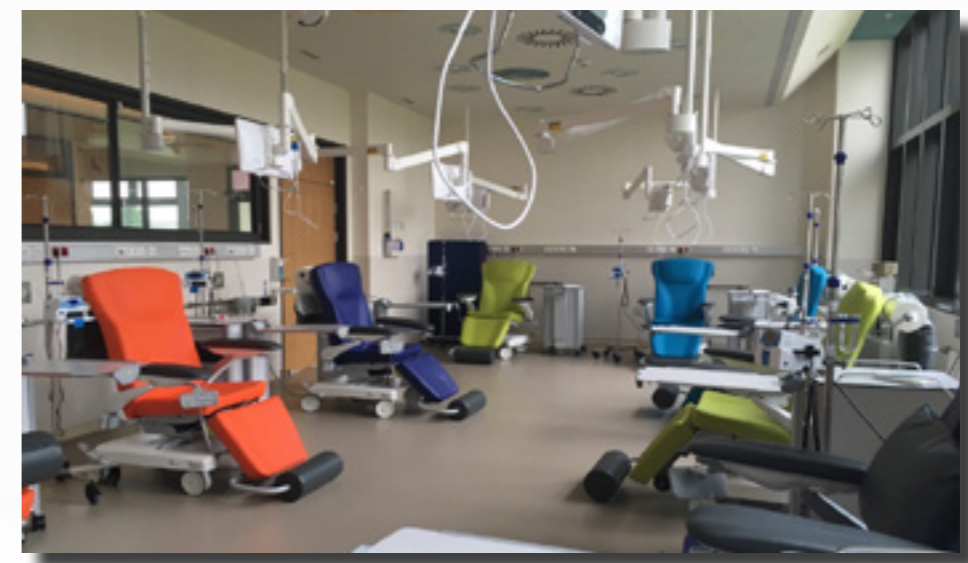

**Comfortable**

**Adaptable** Various options available

**Easy to use** variable electric height

**Versatile** Adapted to different units



## VARIABLE-HEIGHT CARE CHAIR CAREXIA FPV

CAREXIA is a versatile treatment chair.

Its soft and welcoming shapes bring comfort, relaxation and serenity to the patient.

Pure and harmonious design is the vector of a positive image for your unit.

CAREXIA is a chair designed for post-surgery rest, blood sampling, care or examinations, chemotherapy or hemodialysis.




# CAREXIA FPU

## DIMENSIONS CAREXIA FPV

	LENGTH	WIDTH	HEIGHT
OVERALL	76.4 cm (folded) 170 cm (open)	73.2 cm 90 cm (with FPOAC3)	130 cm
SEAT	50 cm	57 cm	55 cm on seat
INCLINE FOLDER	85 to 60°		
VARIABLE HEIGHT			54 to 78 cm

SAFE WORKING LOAD: 200 kg

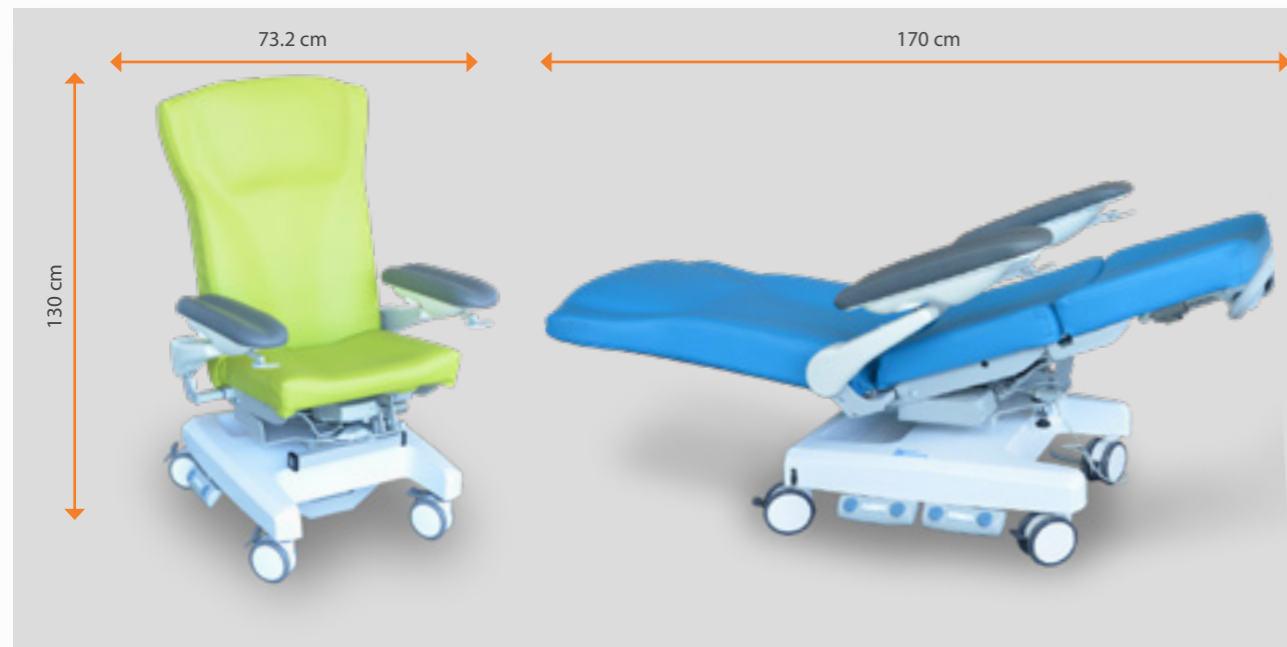

NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

### Height variable electrical care chair

- Heavy duty design. Steel frame with powder coating ABS cover
- Electric variable height from 54 to 78 cm
- Ergonomic shaped backrest and seat. PVC coated jersey
- High resilience foam 35 kg/m3
- Integrated headrest
- Backrest adjustable from 85 to 60° with gas spring assistance (85 to 6° with legrest option)
- Trend position with optional legrest
- 4 self lockable casters with double running wheels
- ABS basket
- Polyurethan armrests, can be folded up



4 self lockable casters with double running wheels



ABS basket



Integrated headrest



Polyurethan armrests, can be folded up


### Faux leather covers : URBAN TEXTILE



ABSINTHE



TOMATO

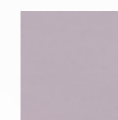

PANSY



FUCHSIA




PLUM



DEEP OCEAN



PETROL

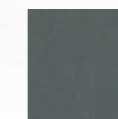

GREY


URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %  
Thickness: 1,2 mm +/- 0,1 mm  
Effective against bacteria (MRSA)  
Phtalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof



FIRE PROOF  
M1

# ACCESSORIES

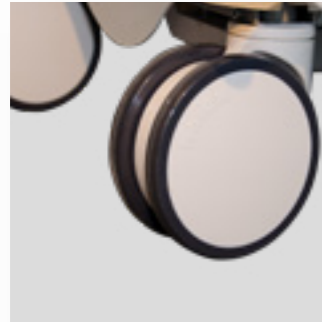

4 castors Ø 125 mm with central lock  
FPORS1



4 double runing wheel castors Ø 150 mm with central lock  
FPORS2



Push bar & Paper roll holder  
FPOBP & FPOSR




Short legrest synchronized with the backrest. Trend position.  
FPORJ



2 holders for IV pole (necessary for TSR2)  
FPOSA



Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)  
TSR2

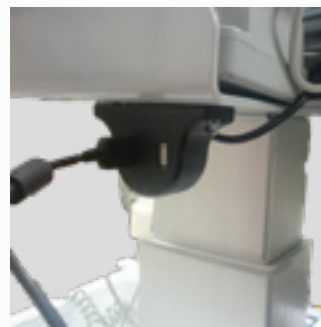

USB Port  
SUP-USB



Comfort pillow (grey)  
ACC-CTG



Armrest comfort (pair)  
FPOAC4



Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair) FPOAC3

REFERENCE	DESCRIPTION
FPORS1	4 castors Ø 125 mm with central lock
FPORS2	4 double runing wheel castors Ø 150 mm with central lock
FPOBP	Push bar
FPORJ	Short legrest synchronized with the backrest. Trend position.
FPOSA	2 holders for IV pole (necessary for TSR2)
TSR2	Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)
TSR4	Removable IV pole, 2-stage, 4 hooks (FPOSA mandatory)
ACC-CTG	Comfort pillow (grey)
FPOAC4	Armrest comfort (pair)
FPOAC3	Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)
FPOSR	Paper roll holder
FPOGC	High comfort upholstery with 2 cm thickness memory foam
FPOBA	Battery 2 Ah
SUP-USB	USB Port

Comfortable

All-electric care chair

Adaptable Various options available

Easy to use Variable height & electric movements

Versatile Adapted to different units



ALL ELECTRIC CARE CHAIR CAREXIA FPVE

CAREXIA is a versatile treatment chair.

Its soft and welcoming shapes bring comfort, relaxation and serenity to the patient.

Pure and harmonious design is the vector of a positive image for your unit.

CAREXIA is a chair designed for post-surgery rest, blood sampling, care or examinations, chemotherapy or hemodialysis.




# CAREXIA FPVE

## CAREXIA FPVE DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	76.4 cm (folded) 177 cm (open) + long leg rest	73.2 cm 90 cm (with FPOAC3)	130 cm
SEAT	50 cm	57 cm	55 cm on seat
INCLINE FOLDER	80 to 6°		
VARIABLE HEIGHT			54 to 78 cm

SAFE WORKING LOAD: 200 kg

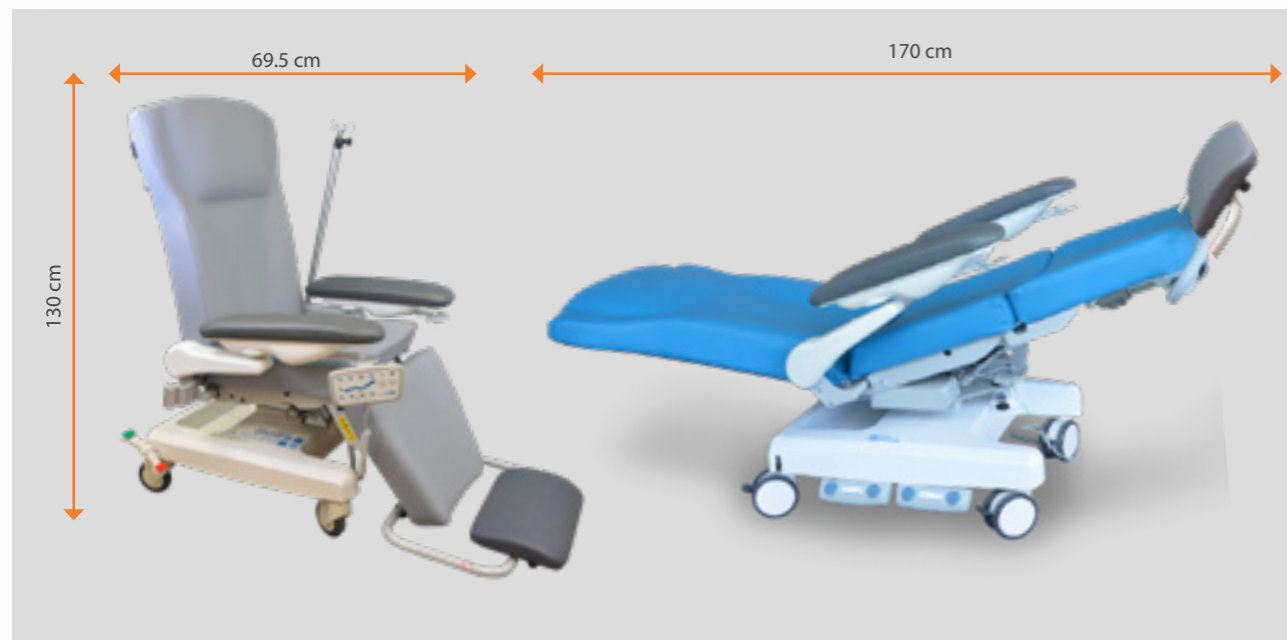

NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

### Height variable electrical care chair

- Heavy duty design. Steel frame with powder coating. ABS cover
- Electric variable height from 54 to 78 cm
- Ergonomic shaped backrest and seat. PVC coated jersey
- High resilience foam 35 kg/m3..
- Integrated headrest
- Electric backrest and legrest
- Trend position
- 4 self lockable casters with double running wheels
- ABS basket
- Polyurethan armrests, can be folded up



4 self lockable casters with double running wheels




Integrated headrest



ABS basket



Trend position

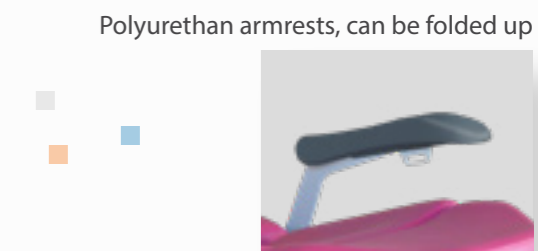

Polyurethan armrests, can be folded up

### Faux leather covers : URBAN TEXTILE



ABSINTHE



TOMATO



PANSY



FUCHSIA




PLUM



DEEP OCEAN



PETROL

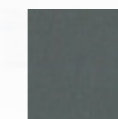

GREY

URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %  
Thickness: 1,2 mm +/- 0,1 mm  
Effective against bacteria (MRSA)  
Phthalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof




FIRE PROOF  
M1



## ACCESSORIES

RÉFÉRENCE	DÉSIGNATION
FPORS1	4 castors Ø 125 mm with central lock
FPORS2	4 double runing wheel castors Ø 150 mm with central lock
FPOBP	Push bar
FPORJL	Long legrest
FPORP	Electric adjustable footrest (only for long legrest) (incompatible with FPOAC4)
FPOCP1	Foot pedal control for relax and trend position
FPOCP2	Foot pedal control for variable height on both sides
FPOSA	2 holders for IV pole (necessary for TSR2)
TSR2	Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)
TSR4	Removable IV pole, 2-stage, 4 hooks (FPOSA mandatory)
ACC-CTG	Comfort pillow (grey)
FPOAC4	Armrest comfort (pair)
FPOAC3	Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)
FPOSR	Paper roll holder
FPOGC	High comfort upholstery with 2 cm thickness memory foam
FPOLAMP	Reading lamp DEL
FPOBA	Battery 2 Ah
SUP-USB	USB Port
PROTEC RJL	Removable Crystal protection for footrest
PROTEC RP	Removable Crystal protection for footrest
SUPTAB	Fixed I-Pad support with antitheft device

# ACCESSORIES

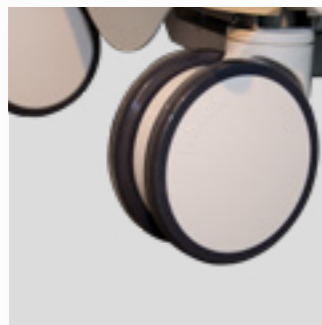

4 castors Ø 125 mm with central lock  
FPORS1



4 double running wheel castors Ø 150 mm with central lock  
FPORS2



Push bar & Paper roll holder  
FPOBP & FPOSR




Long legrest & Electric adjustable footrest (only for long legrest)  
FPORJL & FPORP



FPOCP1 & FPOCP2




2 holders for IV pole (necessary for TSR2)  
FPOSA

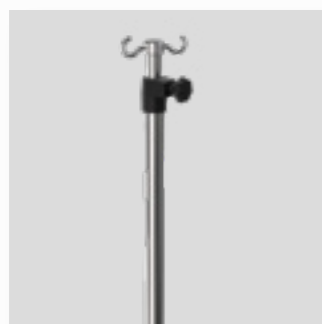

Removable IV pole, 2-stage, 2 hooks (FPOSA mandatory)  
TSR2



Comfort pillow (grey)  
ACC-CTG

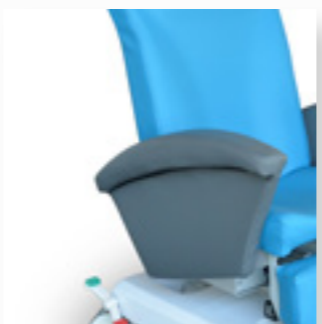

Armrest comfort (pair)  
FPOAC4




Anatomic shaped armrest adjustable in height (2 positions) and in all directions. Retractable behind the backrest for lateral transfer of the patient. (pair)  
FPOAC3



Reading lamp DEL  
FPOLAMP



USB Port  
SUP-USB

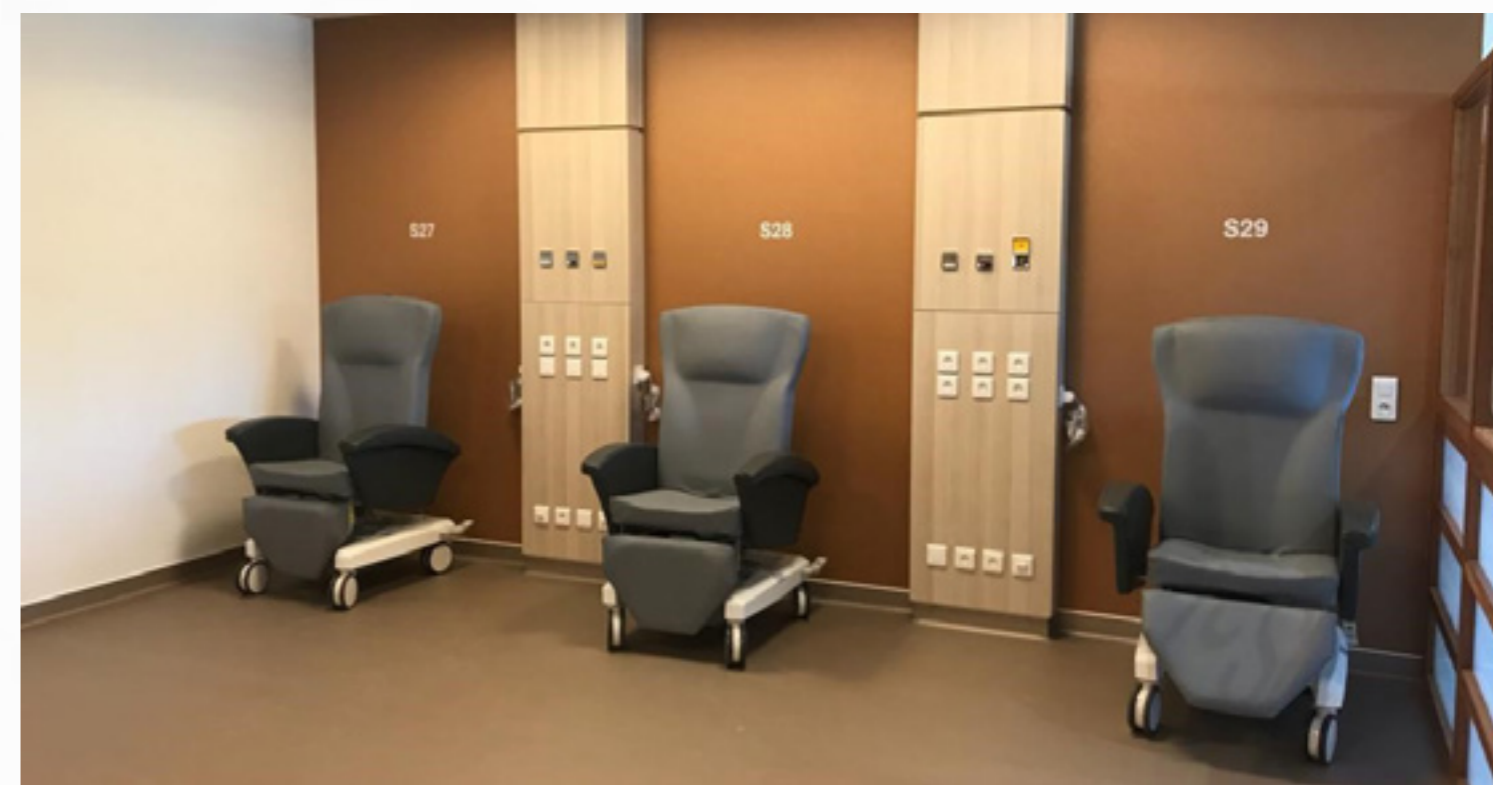





# TIM MANUAL RECLINER

**Versatile:** Adapts to different hospital services

**Easy to use** Manually programmed positions

**Adaptable** Various options available

**Customizable** Choose your wheelbases and armrests

The TIM is a manual recliner that can be adapted to different hospital services.

Its ergonomic design and hotel style provide the patient with exceptional comfort.

The TIM is a versatile recliner with all of the functions required for procedures with medical equipment.

The recliner has been designed for routine use, as well as to simplify regular maintenance.




# T I M

## TIM DIMENSIONS

	LENGTH	WIDTH	HEIGHT	DEPTH
OVERALL	180 cm	72 cm 80 cm (with aluminum armrests)	118 cm	88 cm
SEAT		55 cm	45 cm	50 cm
LEG REST	40 cm	55 cm		
BACKREST	80 cm	58 cm		
HEADREST		55 cm		

SAFE WORKING LOAD: 130 kg

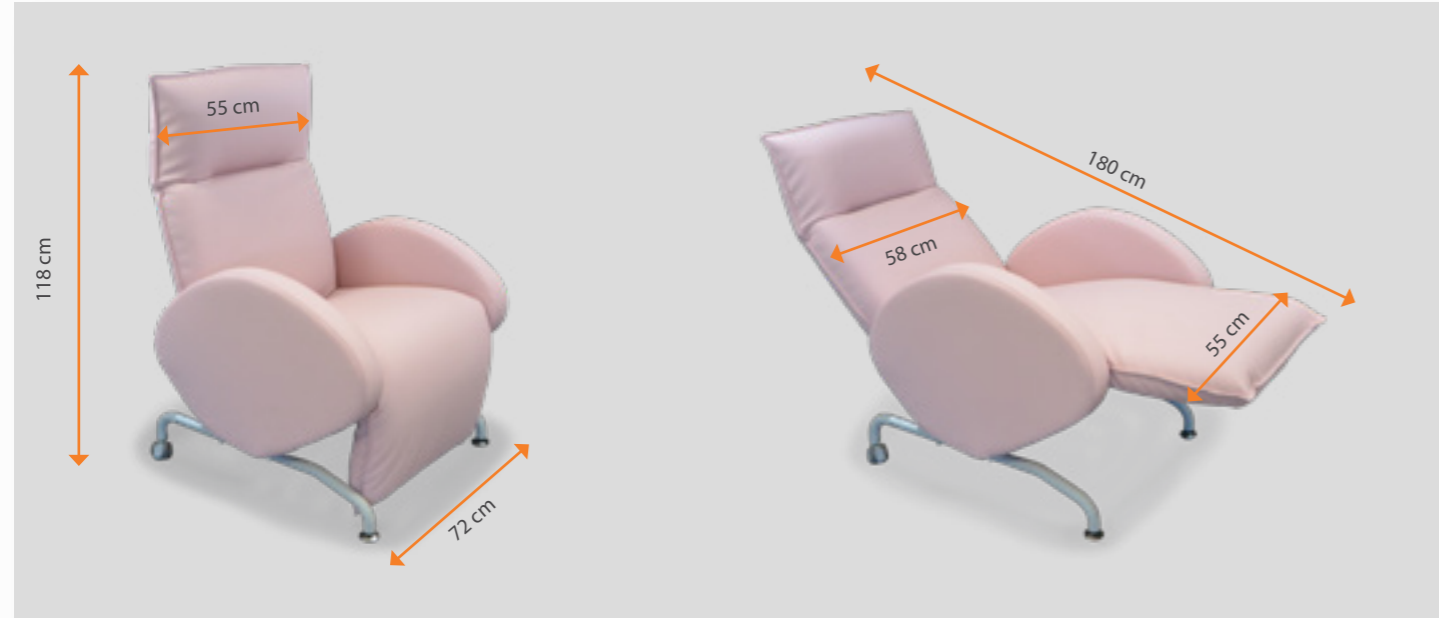

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Day Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.


Manual recliner with steel frame. All cushions are retractable for cleaning or repairs.

- Backrest and leg rest synchronize until a comfortable position is reached.
- The backrest reclines independently to the recovery position.
- Reinforced steel tubular frame
- Adjustable headrest using a rack and pinion system
- Seat on helical springs
- Seat cushions and backrest are made from 35 kg density HR foam
- Detachable comfortable armrests
- M2 flame retardant coating
- Frontal seat
- "Leg" section under the seat opens fully
- Tubular frame: Front: 2 adjustable pads, Rear: 2 wheels with brakes (Ø50 mm)



Tubular frame: Front: 2 adjustable pads, Rear: 2 wheels with brakes (Ø50 mm)

"Leg" section under the seat opens fully

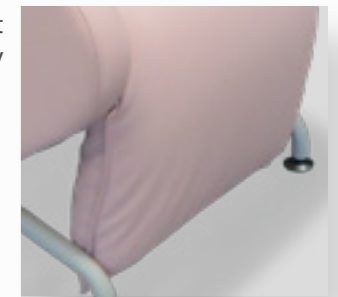

Synchronized backrest and leg rest

Detachable comfortable armrests



# T I M

Adjustable headrest using a rack and pinion system



Seat cushions and backrest are made from 35 kg density HR foam



M2 flame retardant coating

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.




Faux leather covers: VALENCIA M2


## SOOTHING



## REFRESHING

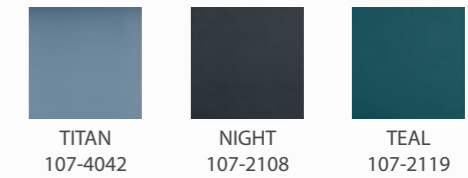

## ENERGIZING

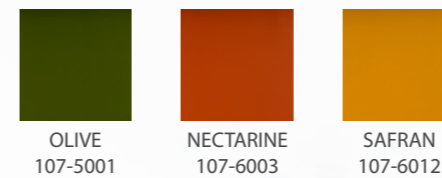

Composition : Vinyl / Urethane Topcoat  
 Thickness : 1,1 mm  
 >300.000 Martindale cycles  
 Anti-stain finish  
 Resistant to sulfide stains  
 Effective against bacteria  
 Phthalate free  
 Anti-static finish



# TIM



## ACCESSORIES

REFERENCE	NAME
FA007	Reading light + USB port
FA021	Fixed aluminum armrests
FA025	Tubular wheelbase with 4 wheels with independent brakes Ø50 mm
FA030	Star wheelbase: 3 wheels with brakes + 2 fixed pads at the front
FA019	Stainless steel base Ø 690 mm
FA031	Cervical roll cylindrical neck pillow
FA034	Legrest crystal protection



# ACCESSORIES




Reading light + USB port  
FA007



Fixed aluminum armrests  
FA021



Cervical roll cylindrical neck  
pillow  
FA031

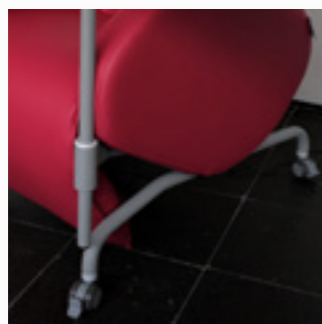

Tubular wheelbase with 4  
wheels with independent  
brakes Ø50mm  
FA025



Stainless steel base  
Ø 690 mm  
FA019



Tubular wheelbase with 2  
rear wheels + 2 front pads  
FA030



Legrest crystal protection  
FA034





# VELVET ELECTRIC RECLINER

**Versatile:** Adapts to different hospital services

**Easy to use** Electrically programmed positions

**Adaptable** Various options available

The VELVET is an electric recliner that can be adapted to different hospital services.

Its ergonomic design and hotel style provide the patient with exceptional comfort.

The VELVET is a versatile recliner with all of the functions required for procedures with medical equipment.

The recliner has been designed for routine use, as well as to simplify regular maintenance.



# VELVET

## VELVET DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	100 cm (folded) 185 cm (open)	82 cm (with ALU armrests) 74 cm (with comfort armrests)	125 cm in seated position
SEAT	50 cm	58 cm	50 cm
LEG REST	40 cm	55 cm	
BACKREST	80 cm	58 cm	
HEADREST		55 cm	

SAFE WORKING LOAD: 160 kg



Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Electric recliner with steel frame. All cushions are retractable for cleaning or repairs.

- Independent electric backrest and leg rest.
- Reinforced steel tubular frame with epoxy coating
- Adjustable headrest using a rack and pinion system
- Seat on combined movement helical springs
- Seat cushions and backrest are made from 35 kg density HR foam
- Armrests with 7 cm of HR foam filling and attached to aluminum frames
- M2 flame retardant coating
- Frontal seat
- "Leg" section under the seat opens fully
- Fixed height
- Star base on wheels ø 75 mm with independent brakes
- Remote control holder




Star base on wheels ø 75 mm with independent brakes

"Leg" section under the seat opens fully



Aluminum armrests with strap for remote control

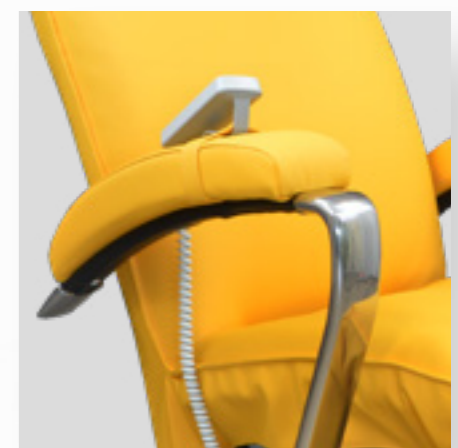

Independent electric backrest and leg rest.

# VELVET

Adjustable headrest using a rack and pinion system



Seat cushions and backrest are made from 35 kg density HR foam



M2 flame retardant coating

Faux leather covers: VALENCIA M2

## SOOTHING




GRAPHITE  
107-4003



TAUPE  
107-0034



SISAL  
107-1048

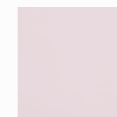

ROSE  
107-2114

## REFRESHING



TITAN  
107-4042



NIGHT  
107-2108



TEAL  
107-2119

## ENERGIZING



OLIVE  
107-5001



NECTARINE  
107-6003



SAFRAN  
107-6012

Composition : Vinyl / Urethane Topcoat  
Thickness : 1,1 mm  
>300.000 Martindale cycles  
Anti-stain finish  
Resistant to sulfide stains  
Effective against bacteria  
Phthalate free  
Anti-static finish



FIRE RATING  
M2

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.



# VELVET



## ACCESSORIES

REFERENCE	DESCRIPTION
FA006	Emergency electric recline function on remote control
FA007	Reading light + USB port
FA008VE	Bilateral holder for IV pole
TSR2	Removable IV pole, 2-stage, 2 hooks (FA008VE mandatory)
TSR4	Removable IV pole, 2-stage, 4 hooks (FA008VE mandatory)
FA011	Examination sheet roll holder
FA015B	Retractable aluminum armrests with 2 positions
FA022	DETACHABLE comfortable "full" armrests (pair) incompatible with FA018V
FA009P	Small Shelf : 42 x 32 cm (FA018V-D or FA018V-G necessary)
FA009G	Large Shelf : 81 x 41 cm (FA018V-D or FA018V-G necessary)
FA018V-D	Support for meal tray on the right side of the patient (compatible with the FA018V-G)
FA018V-G	Support for meal tray on the right side of the patient (compatible with the FA018V-D)
FA025	Tubular wheelbase with 4 wheels with independent brakes Ø50 mm
FA033	Tubular wheelbase with 2 rear wheels + 2 front pads
FA019	Stainless steel base Ø 690 mm
FA031	Cervical roll cylindrical neck pillow
FA034	Legrest crystal protection

# ACCESSORIES




Reading light + USB port  
FA007




Cervical roll cylindrical neck  
pillow  
FA031



Emergency electric  
recline  
on remote control  
FA006

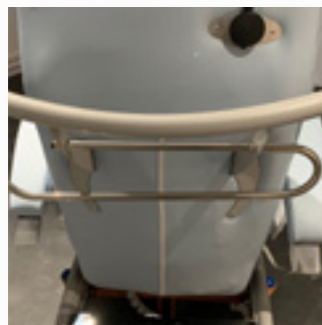


Examination sheet roll holder  
FA011

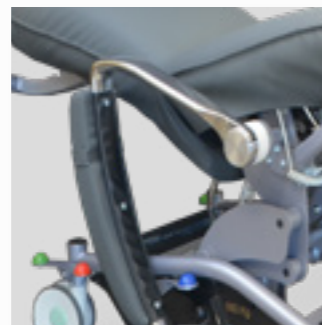


Retractable aluminum  
armrests with 2 positions  
FA015B



DETACHABLE comfortable  
"full" armrests (pair) incom-  
patible with FA018V  
FA022

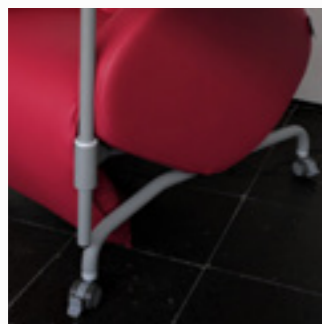

Tubular wheelbase with 4  
wheels with independent  
brakes Ø50mm  
FA025

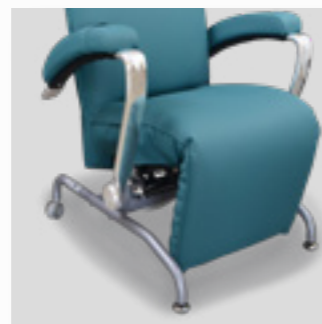


Tubular wheelbase with 2  
rear wheels + 2 front pads  
FA033



Stainless steel base  
Ø 690 mm  
FA019



Small Shelf : 42 x 32 cm  
(FA018V necessary)  
FA009P

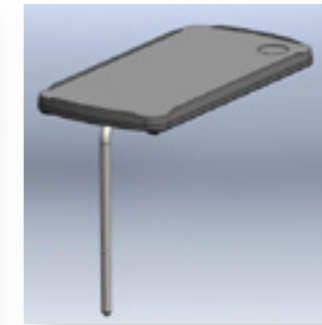

Large Shelf : 81 x 41 cm  
(FA018V necessary)  
FA009G



Velvet Shelf support  
FA018V



Bilateral holder for IV pole  
FA008VE



Legrest crystal protection  
FA034





SWANN

ELECTRIC RECLINER

**Versatile** Adapts to different hospital services

**Easy to use** Electrically programmed positions

**Adaptable** Various options available



The SWANN is an electric healthcare recliner that can be adapted to different hospital services.

Its ergonomic design and hotel style provide the patient with exceptional comfort.

The SWANN is a versatile recliner with all of the functions required for procedures with medical equipment.


The recliner has been designed for routine use, as well as to simplify regular maintenance.

# SWANN

## SWANN DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	100 cm (folded) 185 cm (open)	80 cm	126 cm in seated position
SEAT	50 cm	58 cm	51 to 78 cm
LEG REST	40 cm	55 cm	
BACKREST	80 cm	58 cm	
HEADREST		55 cm	

SAFE WORKING LOAD: 160 kg

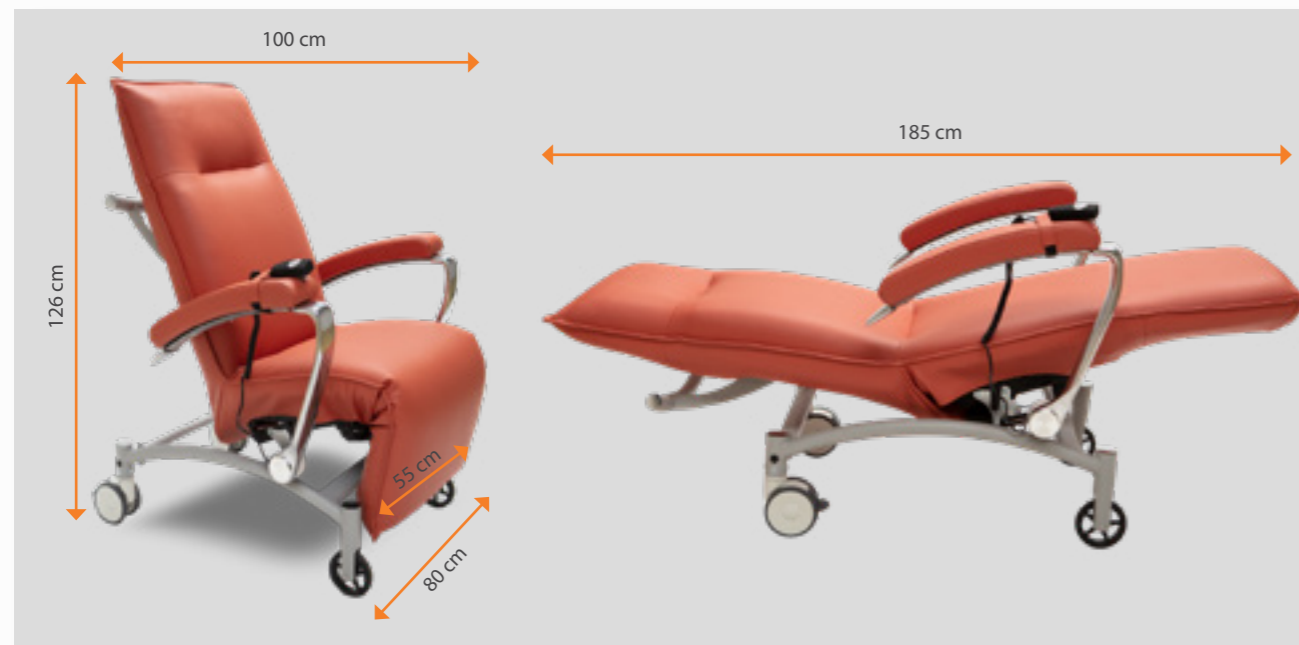

Our range of REST AND HEALTHCARE RECLINERS is adapted for the Day Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Electric recliner with steel frame. All cushions are retractable for cleaning or repairs.

- Fixed height tubular frame
- 2 swiveling rear wheels Ø 125 mm with independent brakes
- 2 fixed front wheels Ø 125 mm
- Remote-controlled electric emergency recline
- PVC-coated jersey.
- M2 flame retardant coating
- Aluminum armrests
- Adjustable headrest using a rack and pinion system
- 2 holders for telescopic IV poles
- Electrically-operated backrest and leg rest.
- Push handle




2 swiveling rear wheels Ø 125 mm with independent brakes

2 fixed front wheels Ø 125 mm



Aluminum armrests with strap for remote control

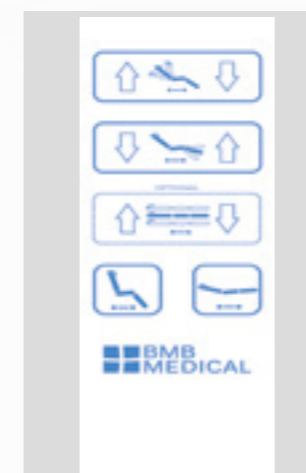

Remote-controlled electric emergency recline

# SWANN



Adjustable headrest using a rack and pinion system

Electrically operated backrest and leg rest.



Push handle



M2 flame retardant coating



2 holders for telescopic IV poles



Our range of REST AND HEALTHCARE RECLINERS is adapted for the Outpatient Surgery, Day Care, Cancer Centre, General Medicine, Oncology and Dialysis units.


Combining comfort, ergonomics and aesthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Faux leather covers: VALENCIA M2


## SOOTHING



## REFRESHING

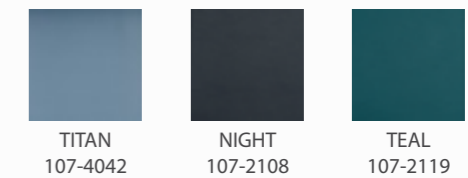

## ENERGIZING

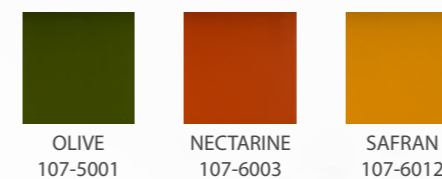

Composition : Vinyl / Urethane Topcoat  
 Thickness : 1,1 mm  
 >300.000 Martindale cycles  
 Anti-stain finish  
 Resistant to sulfide stains  
 Effective against bacteria  
 Phthalate free  
 Anti-static finish





## ACCESSORIES

REFERENCE	DESCRIPTION
FA007	Reading light + USB port
FA008	IV pole
FA010	Battery
FA011	Examination sheet roll holder
FA012	Anatomic shaped armrest - Dimensions 400/200/40 mm (FA012S + FA015 or FA015B required)
FA012S	Support for Anatomic shaped armrest - The pair (FA015 or FA015B required)
FA014	Electrically-powered adjustable height
FA014R	Variable height with lift
FA017	Bilateral foot pedal for reclining (Electric height adjustment necessary FA014)
FA015	Retractable armrests (pair) with 3 positions (FA014 obligatory)
FA015B	Retractable aluminum armrests with 2 positions
FA009P	Small Shelf : 42 x 32 cm (FA018S-D or FA018S-G necessary)
FA009G	Large Shelf : 81 x 41 cm (FA018S-D or FA018S-G necessary)
FA018S-D	Support for meal tray on the right side of the patient (compatible with the FA018S-G)
FA018S-G	Support for meal tray on the right side of the patient (compatible with the FA018S-D)
FA026	Centralized rear wheel locking via a pedal Ø 150 mm with double loading rollers
FA029	2 swiveling double loading rollers to the front with brakes Ø 100 mm
FA027	Floor-facing light beams
FA031	Cervical roll cylindrical neck pillow
FA034	Removable legrest crystal protection

# ACCESSORIES




Reading light + USB port  
FA007



Cervical roll cylindrical neck  
pillow  
FA031



Battery  
FA010

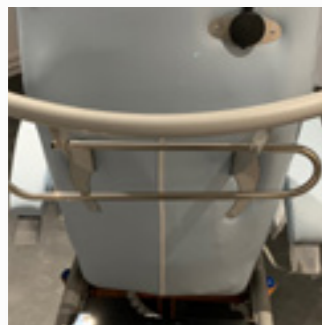


Examination sheet roll holder  
FA011



Phlebotomy armrest  
+ Stand  
FA012 - FA012S



Electrically-powered adjust-  
able height  
FA014

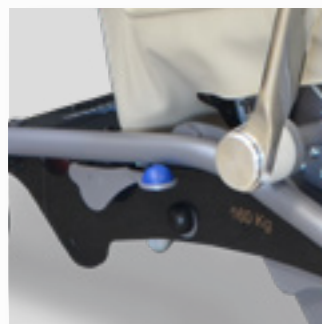

Bilateral foot pedal  
for reclining  
FA017




Retractable armrests (pair)  
with 3 positions  
(FA014 obligatory)  
FA015



Retractable aluminum  
armrests with 2 positions  
FA015B



Small Shelf : 42 x 32 cm  
(FA018S necessary)  
FA009P

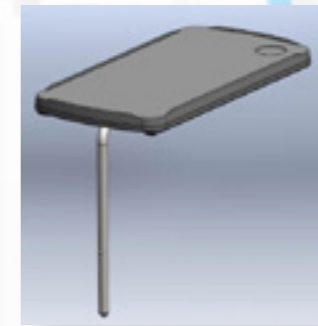

Large Shelf : 81 x 41 cm  
(FA018S necessary)  
FA009G



SwannShelf support  
FA018S



2 swiveling double load-  
ing rollers to the front  
with brakes Ø 100 mm  
FA029



Removable legrest crystal protection  
FA034



Centralized rear wheel lock-  
ing via a pedal Ø 150 mm  
with double loading rollers  
FA026




# URBAN TEXTILE M1


## CAREXIA RANGE




ABSINTHE

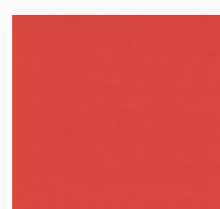


TOMATO

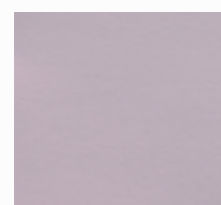


PANSY

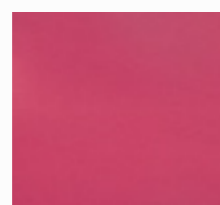

FUCHSIA

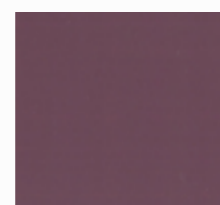


PLUM

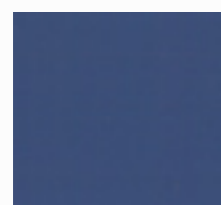

DEEP OCEAN



PETROL

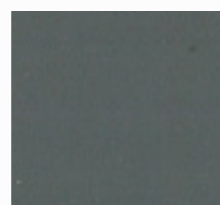

GREY



## COLOR CHART



### CLASSIFICATION


M1

Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %  
Thickness: 1,2 mm +/- 0,1 mm



Effective  
against bacteria  
(MRSA)

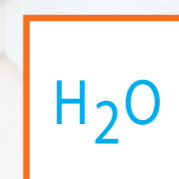

Solvent-free



Phtalate-free



Martindale  
abrasion 100000  
cycles




Waterproof



# VALENCIA M2


## SOOTHING

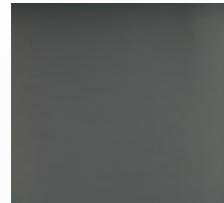


GRAPHITE  
107-4003



TAUPE  
107-0034

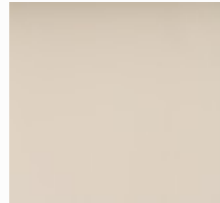

SISAL  
107-1048

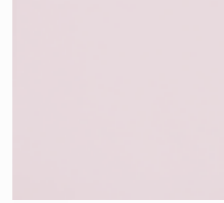


ROSE  
107-2114

## REFRESHING



TITAN  
107-4042

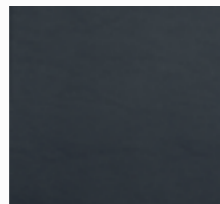

NIGHT  
107-2108



TEAL  
107-2119

## ENERGIZING



OLIVE  
107-5001



NECTARINE  
107-6003



SAFRAN  
107-6012

## COLOR CHART



FIRE RATING  
M2

COMPOSITION : Vinyl / Urethane Topcoat  
THICKNESS : 1,1mm

### FINISHINGS AND TREATMENTS



Effective  
against bacteria  
(MRSA)



Antistatic  
finish



Phtalates free



>300000  
Martindale cycles



Stain and sulfide  
resistant finish

A medical professional must always rely on his or her own professional clinical judgment when deciding whether to use a product(s) when treating a particular patient.  
A healthcare professional anyway, refers to the package insert, product label and/or instructions for use before using any BMB MEDICAL product.  
The products depicted are CE marked in accordance with EU Regulations and Directives.  
The products are Class 1 Medical Devices.



Tél : +33 (0)4 74 08 71 71

Fax : +33 (0)4 74 08 71 72

[info@bmb-medical.com](mailto:info@bmb-medical.com)

100 Av. du Formans - CS 11001

01604 TREVoux Cedex

FRANCE

[www.bmb-medical.com](http://www.bmb-medical.com)