


EMERGENCY AND REST CHAIR KYELLE



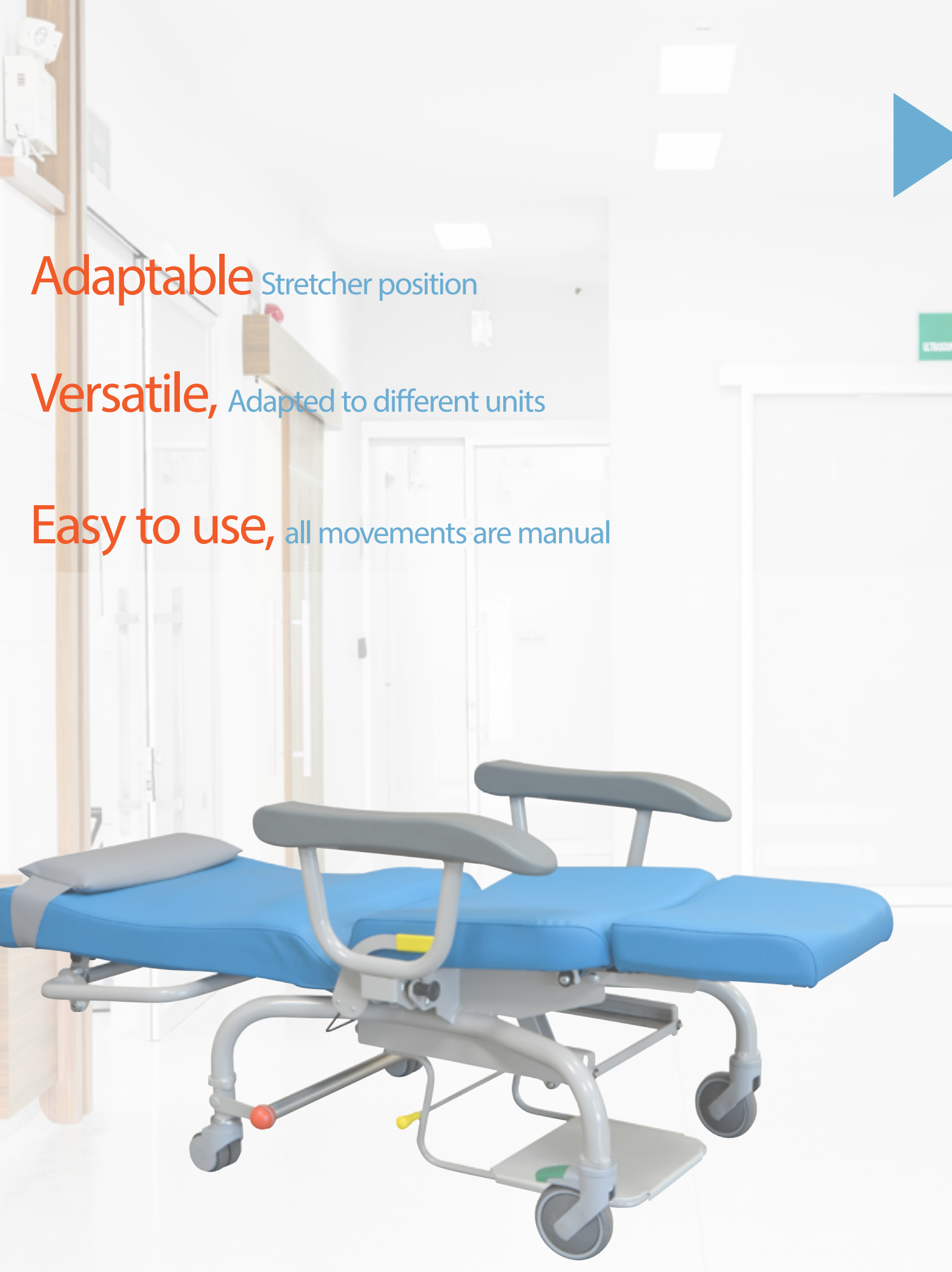

Adaptable Stretcher position

Versatile, Adapted to different units

Easy to use, all movements are manual



M A N U A L C H A I R K Y E L L E

The KYELLE is a manual chair, for reception and transfer, suitable for emergency units and other hospital departments.

Its robust and ergonomic design makes it a real multi-purpose product.

The KYELLE is an emergency reception chair but also a rest chair.




KYELLE

KYELLE DIMENSIONS

| | LENGHT | WIDTH | HEIGHT |
|---------|---|-------|--------|
| OVERALL | 76 cm in sitting position 161 cm in lying position | 73 cm | 120 cm |
| SITING | | 52 cm | 48 cm |

SAFE WORKING LOAD : 200 Kg

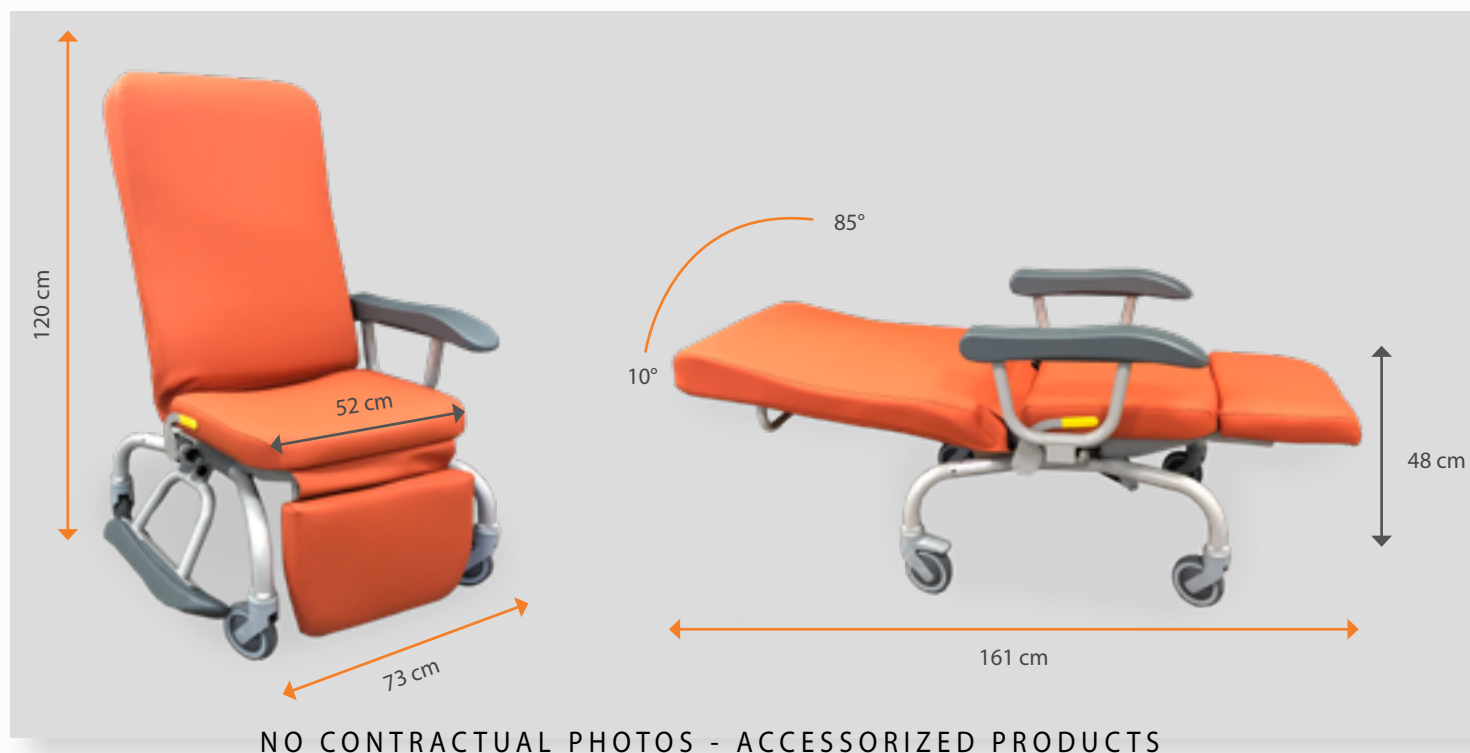

NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Manual chair with epoxy coated steel structure

- Epoxy coated steel frame.
- Assembled on 4 pads (2 back pads with integrated casters)
- Unclippable seat. Ergonomically shaped backrest. M1 fireproof PVC covering.
- Range of colors URBAN.
- Foam upholstery thickness 8cm
- Fixed backrest inclined at 85°.
- Fixed polyurethane armrests
- Weight: 27 to 34 kg depending on model
- Safe working load: 200 kg.
- Maximum patient weight: 135 kg.
- Warranty: 2 years



Fixed polyurethane armrests



Foam upholstery thickness 8cm

Faux leather covers : URBAN M1



URBAN
Fine grain, mat and gentle to the touch.
Rot-proof cotton base.


Composition : PVC 85 % / Coton 15 %
Thickness: 1,2 mm +/- 0,1 mm
Effective against bacteria (MRSA)
Phtalate-free
Solvent-free
Martindale abrasion 100 000 cycles
Waterproof



ACCESSORIES



4 casters of 125mm diameter:
2 free casters at the front, 2
rear casters with independent
brakes
KYRF

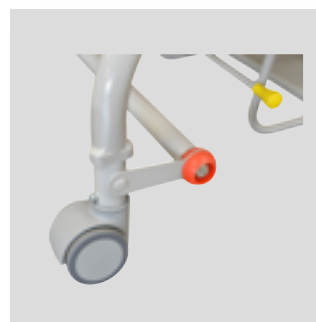

2 free casters at the front of
125mm diameter, 2 centrally
braked rear casters of 75mm
diameter
KYRS



Ergonomically shaped
polyurethane retractable
armrests.
KYAC2




Backrest adjustable from 85
to 10° assisted by gas spring.
KYRB



Leg rest synchronized with the
back rest adjustable from 0 to
100°. (KYRB mandatory)
KYRJ1



Independent leg rest
adjustable from 0 to 100°
assisted by
gas spring.
KYRJ2

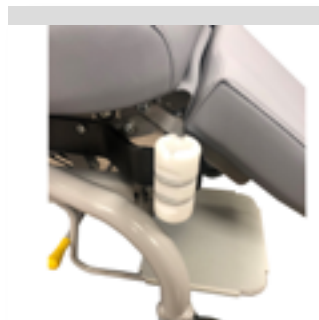

Support for meal tray on the
right/left side of the patient
KYSPR-D or KYSPR-G



FA009P
Small Shelf : 42 x 32 cm
(KYSPR-D or KYSPR-G mandatory)




FA009G
Large Shelf : 81 x 41 cm
(KYSPR-D or KYSPR-G mandatory)



Reversible push bar
(sitting and lying position)
(KYRF or KYRS mandatory)
KYBP2



Fixed push bar
(KYRF or KYRS mandatory)
KYBP1

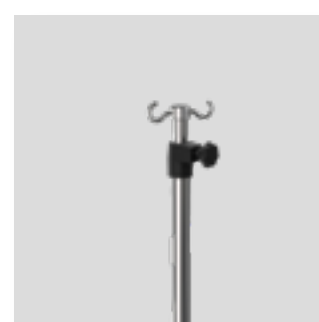

Removable IV Pole
(KYSTR mandatory)
TSR2

| REFERENCE | DESCRIPTION |
|-----------|---|
| KYRF | 4 casters of 125mm diameter : 2 free casters at the front, 2 rear casters with independent brakes |
| KYRS | 2 free casters at the front of 125mm diameter, 2 centrally braked rear casters of 75mm diameter |
| KYRD | Front directional wheel (KYRF or KYRS mandatory) |
| KYAC2 | Ergonomically shaped polyurethane retractable armrests. |
| KYRB | Backrest adjustable from 85 to 10° assisted by gas spring. |
| KYRJ1 | Leg rest synchronized with the back rest adjustable from 0 to 100°. (KYRB mandatory) |
| KYRJ2 | Independent leg rest adjustable from 0 to 100° assisted by gas spring. |
| KYRP | Retractable footrest (KYRF or KYRS mandatory) |
| KYBP1 | Fixed push bar (KYRF or KYRS mandatory) |
| KYBP2 | Reversible push bar (sitting and lying position) - (KYRF or KYRS mandatory) |
| KYCC | Comfort headrest cushion |
| TSR2 | Removable IV pole, 2-stage, 2 hooks (KYSTR mandatory) |
| TSR4 | Removable IV pole, 2-stage, 4 hooks (KYSTR mandatory) |
| KYSTR | IV Pole holder (on both sides of the product) |
| KYSPR-D | Support for meal tray on the right side of the patient (compatible with the KYSPR-G) |
| KYSPR-G | Support for meal tray on the left side of the patient (compatible with the KYSPR-D) |
| FA009P | Small Shelf : 42 x 32 cm (KYSPR-D or KYSPR-G mandatory) |
| FA009G | Large Shelf : 81 x 41 cm (KYSPR-D or KYSPR-G mandatory) |


URBAN TEXTILE M1




ABSINTHE




TOMATO

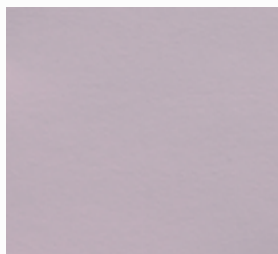


PANSY

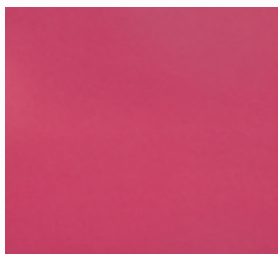

FUCHSIA

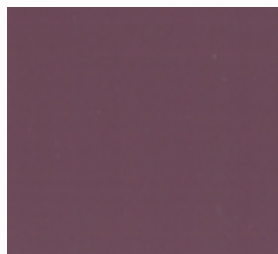


PLUM

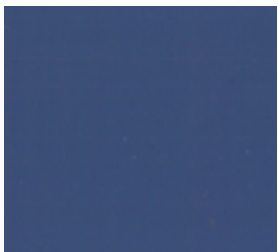

DEEP OCEAN



PETROLE

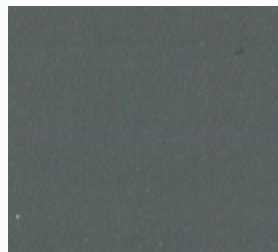

GREY

COLOR CHART



CLASSIFICATION

M1


Fine grain, mat and gentle to the touch.
Rot-proof cotton base.

Composition : PVC 85 % / Coton 15 %

Thickness : 1,2 mm +/- 0,1 mm



Effective
against
bacteria(MRSA)

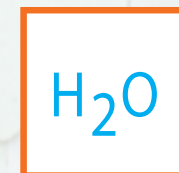

Solvent-free



Phtalate-free



Martindale
a b r a s i o n
1000000 cycles



Waterproof

A medical professional must always rely on his or her own professional clinical judgment when deciding whether to use a product(s) when treating a particular patient.

A healthcare professional anyway, refers to the package insert, product label and/or instructions for use before using any BMB MEDICAL product.

The products depicted are CE marked in accordance with EU Regulations and Directives.

BMB MEDICAL – 100 Avenue du Formans – CS 11001 – 01604 TREVOUX CEDEX – FRANCE

Tel : +33 474 08 71 71 – Fax : +33 474 08 71 72 – email : info@bmb-medical.com

Sarl au capital de 600 000 € - RCS Bourg en Bresse 89B652

SIRET : 950 413 229 00036 – Code NAF : 3250A – Code TVA : FR68 950 413 229

www.bmb-medical.com



Tél : +33 (0)4 74 08 71 71

Fax : +33 (0)4 74 08 71 72

info@bmb-medical.com

100 Av. du Formans - CS 11001
01604 TREVOUX Cedex
FRANCE

www.bmb-medical.com