

ERGOLYS

DIALYSIS BED-CHAIR



Our primary objective is your **satisfaction**. Our **requirements** are as exacting as we would expect as a customer.

By supporting you from defining the product and throughout its life, we guarantee **high-performance** products, and a **responsive** and **high-quality service**.

We try, with humility, to **understand** your **use**, your **constraints** and to redefine the circuit of your patients to **improve** your **working conditions** and the **efficiency** of your **organization**.

We base our relationships on **shared trust**, making it our priority to be by your side for the **long term**.

RELIABILITY

High-quality components,
Long-term use,
A reduced environmental footprint



Our environmental awareness has led us towards an eco-friendly design.

BMB MEDICAL is concerned about the ecological impact of its products at each stage of their life cycle.

This commitment is reflected in actions on packaging, the "4 R" policy (Reduce, Replace, Recycle, Rethink) and in our efforts on materials. Our products are designed and approved in accordance with standard IEC60601-1-9 that defines the requirements for eco-friendly design.



Our columns are maintenance-free. Our sleeping surfaces are made from an extremely resistant and rot-proof material. Our corner protections are resistant to shocks and aging. Choosing BMB MEDICAL is a guarantee of quality all while offering significant savings on your maintenance budget.

FOCUS DIALYSIS BED-CHAIR



The advantageous comfort of a chair, as comfortable as a bed and as functional as a stretcher. Suitable for any clinical situation, it also accommodates large patients.

The project is an innovative furnishing solution, reduces the psychological discomfort of the patient.

The bed-chair can take different positions, thanks to 5 low voltage motors controlled by a push-button panel easily accessible to the patient.

The study of the details and the technological research reveal practical functional solutions. Perfectly suited to accommodate patients of all conformations, it ensures positive emotions and unprecedented relaxation.

It favors any type of therapeutic treatment in perfect postural conditions.

Choose relaxation, choose ERGOLYS



Our range of **DIALYSIS BED-CHAIR** is adapted to dialysis services.

Combining **comfort** and **ergonomics**, our products are designed to satisfy and meet the needs of patients and caregivers.

These products have been developed by BMB Médical in order to fill the "gaps" in **ergonomics**, comfort or current use of furniture and medical equipment.



C O N T E N T S



E R G O L Y S



C O L O R C H A R T

Hybrid, the comfort of a bed and the functionality of a chair
Adaptable, its functions allow the transfer of patients
Comfortable, it favors the postural position

BED-CHAIR **ERGOLYS**

ERGOLYS is an electric BED-CHAIR designed for DIALYSIS services

Its ergonomic design and comfort make it a true hybrid product. It combines the comfort of a bed with the functionalities of a chair.

The ERGOLYS is a multipurpose chair-bed offering all the functions necessary for the well-being and transfer of patients.

Its composition has been thought to facilitate the current use of the chair-bed, but also to simplify its usual maintenance.



ERGOLYS

ERGOLYS DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	163 cm armchair position 225 cm lying position 235 cm (lying position, footrest retracted)	110 cm	62 to 102 cm

PRODUCT WEIGHT : 125 Kg
SAFE WORKING LOAD : 250 Kg



NON CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

ERGOLYS fits perfectly into a dialysis department.

Combining comfort, ergonomics and aesthetics, our products are designed to satisfy and meet the needs of patients and caregivers

Bed-Chair with electric movements

- Electric leg lift 0 to 25°
- Electric slope 20°
- Electric bust lift 0 to 75°
- 3 section bed 70 x 215 cm
- Ultra-comfort mattress: 3 preformed zones: HR40 / PE 40 / Visco memory foam / PU coated
- Chair-Bed with electric movements
- Set of low frame high frame characterized by armrests
- Adjustable headrest by racks
- 5 possible positions: Reception, comfort position, bed, emergency slope, CPR function
- Protective shell of the frame and bed cover
- 4 wheels with independent brakes Ø 125 mm
- Retractable armrests with multidirectional adjustment 55 x 17 cm (VALENCIA colors)
- VALENCIA color range (only for accessories)




Retractable armrests with multidirectional adjustment



4 wheels with independent brakes Ø 125 mm (VALENCIA color)



Trendelenburg 20°

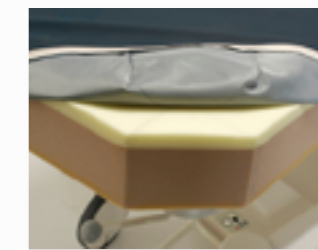

Ultra-comfort mattress



12 boutons dont 5 fonctions automatiques



Electric footrest
ERORPE



Adjustable head cushion
EROCT



Emergency foot switch with
emergency foot release on
both sides
ERODU



Wheels Ø 125 mm with bilateral central-
ized braking, one of which is directional
ERORS1



Flex remote control (satellite)
EROSAT




ACCESSOIRES

RÉFÉRENCE	DÉSIGNATION
ERORPE	Electric footrest 25 cm with retractable foot cushion (VALENCIA color) (EROSAT mandatory)
EROCT	Adjustable head cushion 70 x 45 cm (color VALENCIA)
ERODU	Emergency foot switch with emergency foot release on both sides
ERORS1	Wheels Ø 125 mm with bilateral centralized braking, one of which is directional
EROSAT	Flex remote control (satellite)
EROBA	LI-ion battery

VALENCIA M2

SOOTHING

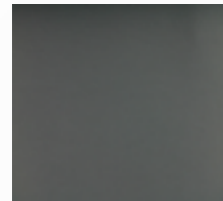


GRAPHITE
107-4003



TAUPE
107-0034



SISAL
107-1048

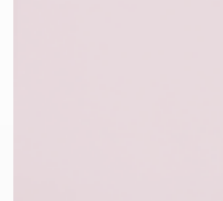


ROSE
107-2114

REFRESHING



TITAN
107-4042

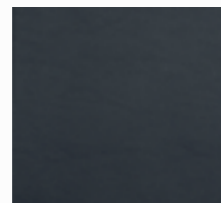

NIGHT
107-2108




TEAL
107-2119

ENERGIZING



OLIVE
107-5001

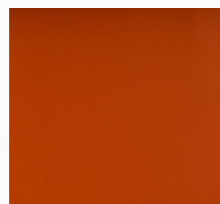

NECTARINE
107-6003



SAFRAN
107-6012

COLOR CHART



FIRE RATING
M2

COMPOSITION : Polyester knitted fabric with fireproof PVC coating on one side
THICKNESS : 1,1mm

FINISHINGS AND TREATMENTS



Effective against
bacteria (MRSA)



Antistatic
finish



Phtalates free



>300000
Martindale cycles



Stain and sulfide
resistant finish

A medical professional must always rely on his or her own professional clinical judgment when deciding whether to use a product(s) when treating a particular patient.

A healthcare professional anyway, refers to the package insert, product label and/or instructions for use before using any BMB MEDICAL product.

The products depicted are CE marked in accordance with EU Regulations and Directives.

The products are Class 1 Medical Devices.



Tél : +33 (0)4 74 08 71 71

Fax : +33 (0)4 74 08 71 72

info@bmb-medical.com

100 Av. du Formans - CS 11001

01604 TREVOUX Cedex

FRANCE

www.bmb-medical.com