

# EXAMINATION COUCH





## CONTENTS


Our range of EXAMINATION COUCH is adapted for the Day Surgery, Day Care, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and esthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

▶ DIVIA

▶ DIVIA MIXTE

▶ COLOR CHART

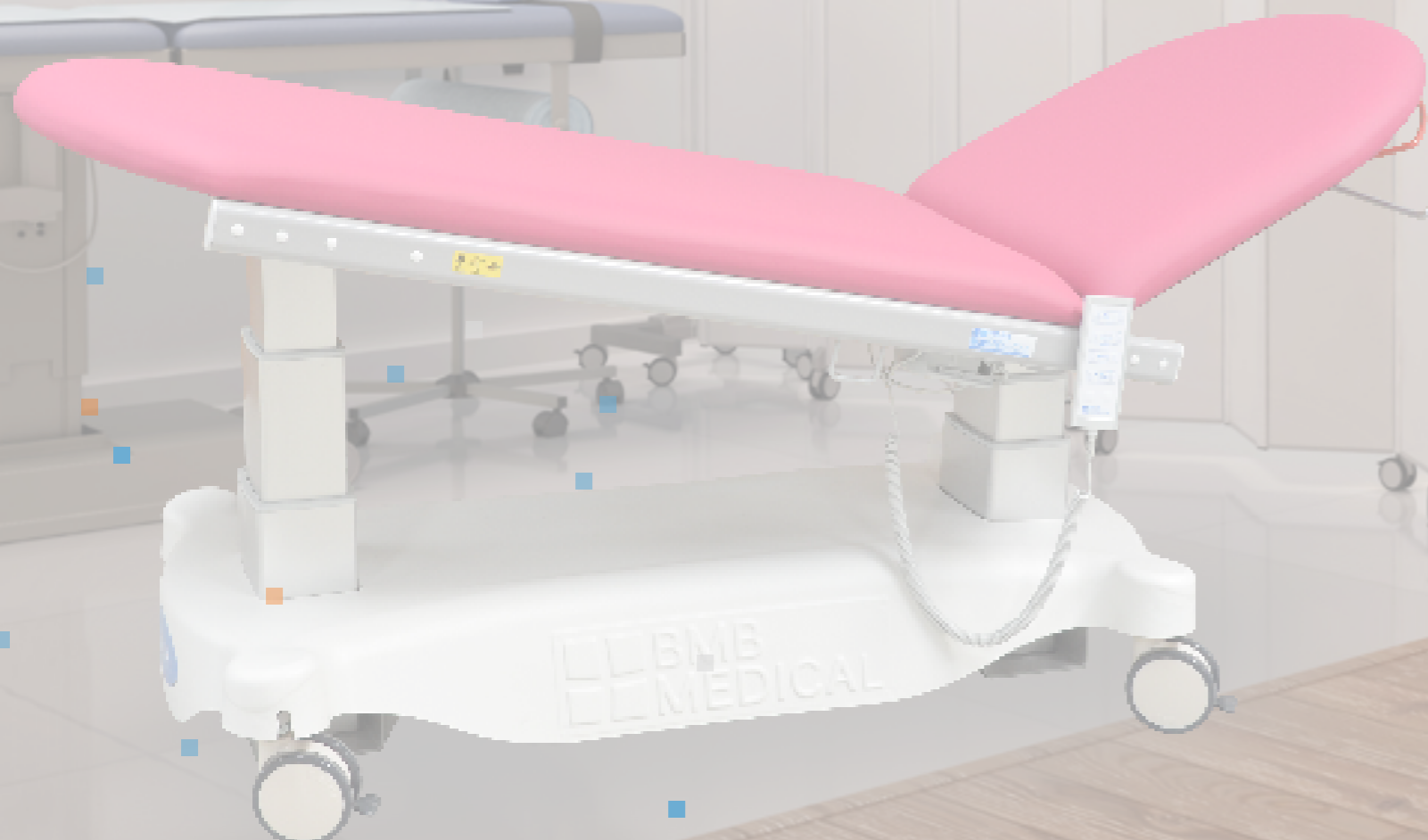

Stability


Intensive use

Easy to use electric height adjustment

## DIVIA VARIABLE-HEIGHT COUCH

DIVIA is an examination couch with electric height adjustment on columns.

The couch is particularly suitable for emergency services and intensive use.

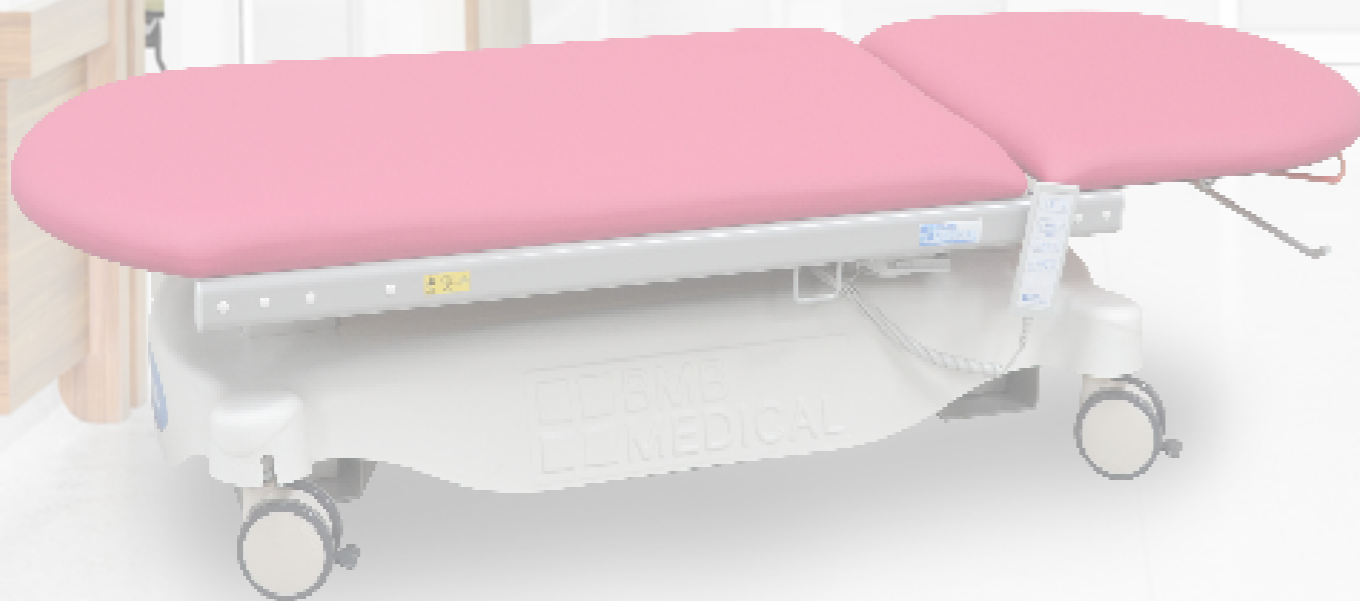

## DIVIA DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	205 cm	73 or 75 cm (depending on model)	47 to 87 cm
SEAT	65 x 200 cm or 75 x 200 cm		
INCLINE FOLDER	0 to 80°		

PRODUCT WEIGHT : 80 or 85 Kg  
SAFE WORKING LOAD: 250 kg



NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of EXAMINATION COUCH is adapted for the Day Surgery, Day Care, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and esthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Examination couch with electric height adjustment on columns.

- Electric height adjustment : 47 to 87 cm
- Electric assisted Trendelenburg / reverse Trendelenburg +/- 15°
- Adjustable headrest 0 - 80° with 2 gas springs
- 2-sectional mattress base : 65x 200 cm or 75 x 200 cm
- Flat upholstery, high density foam 35 kg/m<sup>3</sup> - Thickness 6,5 cm
- Jersey with PVC coating, antibacterial and anti-fire (M1) treatment
- Steel structure, epoxy coated (light grey)
- Frame with 2 electric columns providing great stability
- Frame protection cover
- Remote control
- 4 lockable wheels Ø 125 mm
- Examination towel dispenser head end
- URBAN color range



Remote control



Electric height adjustment : 47 - 87 cm



Examination towel dispenser head end




Electric assisted Trendelenburg / reverse Trendelenburg (+/- 15°)

Faux leather covers : URBAN TEXTILE



URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.  
Composition : Mass-processed PVC 86%, Cotton 14%.  
Thickness : 1,04 mm +/- 10%  
Effective against bacteria (MRSA)  
Phthalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof

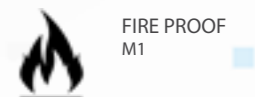

# ACCESSORIES



Antistatic system  
DVRA



Casters Ø 125 mm  
with central brake  
DVRS1



Double running castors Ø  
150 mm with central brake  
DVRS2



Electric adjustable  
backrest  
DVRBE



Wireless remote control  
(in addition of the wired  
remote control)  
DVTSF




Foot pedal control (in addition of the  
wired remote control) If the options  
DVTAP and DVTSF are chosen, the  
wired remote control is removed)  
DVTAP



Drop-down side rails (pair)  
DVBR



Holder for clamps (pair)  
DVSE

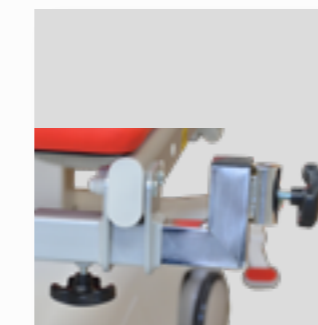

Fixation clamps for heels, armrest  
for blood sampling or leg sup-  
port (pair)  
(DVSE required)  
DVETU



Heels (pair)  
(DVETU required)  
DVET



Leg support (pair)  
(DVETU required)  
DVSJ



Removable 2-stage IV pole, 2 hooks  
(DVTSR required)  
TSR2

REFERENCE	DESCRIPTION
DVRA	Anti-static system
DVRS1	Casters Ø 125 mm with central brake
DVRS2	Double running castors Ø 150 mm with central brake
DVRBE	Electric adjustable backrest
DVTSF	Wireless remote control (available from november)
DVTAP	Foot pedal control (instead of the manual remote control)
DVBR	Drop-down side rails (pair)
DVSE	Holder for clamps (pair)
DVETU	Fixation clamps for heels, armrest for blood sampling or leg support (pair) (DVSE required)
DVET	Heels (pair) (DVETU required)
GT	Anatomic shaped armrest (DVETU necessary)
DVSJ	Leg support (pair) (DVETU required)
DVHSJ	Protection cover for leg rest (pair)
DVTSR	Clamps for IV pole holder (pair) (DVSE necessary)
TSR2	Removable 2-stage IV pole, 2 hooks (DVTSR required)
TSR4	Removable IV pole, 2-stage, 4 hooks (DVTSR required)
ACC-CTG	comfort pillow (grey)
DVPC7	Crystal protection cover for examination couch width 65cm
DVPC8	Crystal protection cover for examination couch width 75cm
DVBBA	Battery
DVLAMP	Examination lamp (DVETU and DVSE required)
EQUI	Equipotential connection
DVRT	Side technical rails (DVSE required)

**Intensive use**

**Versatile** suitable for gynecological examinations

**Easy to use** Electric backrest and height adjustment



## GYNECOLOGICAL COUCH DIVIA MIXTE

DIVIA is a gynecological examination couch with electric height adjustment on columns.

It has a second electrically operated intermediate backrest for gynecological examinations.



# DIVIA MIXTE

## DIVIA MIXTE DIMENSIONS

	LENGTH	WIDTH	HEIGHT
OVERALL	190 cm	73 to 75 cm (depending on model)	47 to 87 cm
SEAT	65 x 178 cm or 75 x 178 cm		
INCLINE FOLDER	0 to 80°		

PRODUCT WEIGHT : 90 Kg  
SAFE WORKING LOAD: 250 kg



NO CONTRACTUAL PHOTOS - ACCESSORIZED PRODUCTS

Our range of EXAMINATION COUCH is adapted for the Day Surgery, Day Care, General Medicine, Oncology and Dialysis units.

Combining comfort, ergonomics and esthetics, our products are designed to accommodate and meet the needs of patients and caregivers.

Examination couch with variable height for gynaecological examinations.

- Electric height adjustment : 47 - 87 cm
- Electric assisted Trendelenburg / reverse Trendelenburg +/- 15°
- Electric adjustable backrest 0 - 80°
- 3-sectional mattress base : 65 x 188 cm or 75 x 188 cm
- Flat upholstery, high density foam 35 kg/m3. Width 6,5 cm
- Steel structure, epoxy coated (light grey)
- Frame with 2 electric columns providing great stability
- Frame protection cover
- Remote control
- 4 lockable double running wheels Ø 125 mm
- Examination towel dispenser head end
- Jersey with PVC coating, antibacterial and anti-fire (M1) treatment
- URBAN color range



Remote control



Electric height adjustment : 46 - 86 cm



Examination towel dispenser head end



Electric assisted Trendelenburg / reverse Trendelenburg (+/- 15°)

### Faux leather covers : URBAN TEXTILE



URBAN  
Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : Mass-processed PVC 86%, Cotton 14%.  
Thickness : 1,04 mm +/- 10%  
Effective against bacteria (MRSA)  
Phtalate-free  
Solvent-free  
Martindale abrasion 100 000 cycles  
Waterproof






Antistatic system  
DVRA




Casters Ø 125 mm  
with central brake  
DVRS1



Double running castors Ø  
150 mm with central brake  
DVRS2

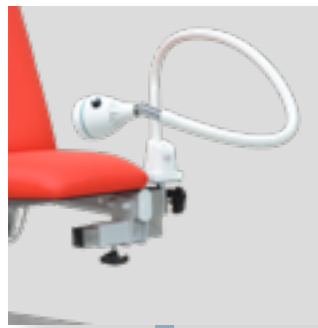

Examination lamp  
(DVETU required)  
DVLAMP

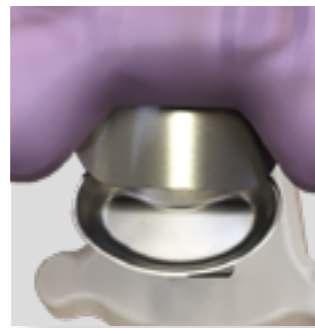

Recovery Tank (grey)  
DVMB



Drop-down side rails (pair)  
DVBR



Holder for clamps (pair)  
DVSE



Fixation clamps for heels, armrest  
for blood sampling or leg sup-  
port (pair)  
(DVSE required)  
DVETU



Heels (pair)  
(DVETU required)  
DVET



Leg support (pair)  
(DVETU required)  
DVSJ



Removable 2-stage IV pole, 2  
hooks  
(DVTSR required)  
TSR2



Comfort pillow (grey)  
ACC-CTG

## ACCESSORIES

REFERENCE	DESCRIPTION
DVRA	Anti-static system
DVRS1	Castors Ø 125 mm with central brake
DVRS2	Double running castors Ø 150 mm with central brake
DVTSF	Wireless remote control (available from november)
DVTAP	Foot pedal control (instead of the manual remote control)
DVBR	Drop-down side rails (pair)
DVSE	Holder for clamps (pair)
DVETU	Fixation clamps for heels, armrest for blood sampling or leg support (pair) (DVSE required)
DVET	Heels (pair) (DVETU required)
DVSJ	leg support (pair) (DVETU required)
DVHSJ	Protection cover for leg rest (pair)
GT	Anatomic shaped armrest (DVETU necessary)
DVTSR	Clamps for IV pole holder (pair) (DVSE necessary)
TSR2	Removable 2-stage IV pole, 2 hooks (DVTSR required)
TSR4	Removable IV pole, 2-stage, 4 hooks (DVTSR required)
ACC-CTG	comfort pillow (grey)
DVBBA	Battery
DVLAMP	Examination lamp (DVETU and DVSE required)
DVMB	Recovery Tank - Can hold up to 0,4 liter
EQUI	Equipotential connection
DVRT	Side technical rails (DVSE required)




# URBAN TEXTILE M1




ABSINTHE




TOMATO

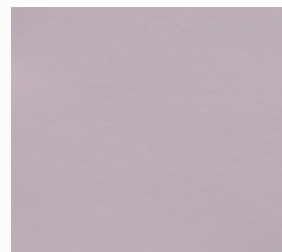


PANSY

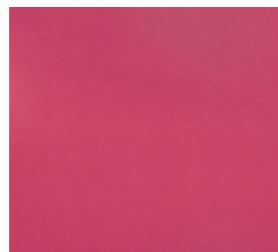

FUCHSIA

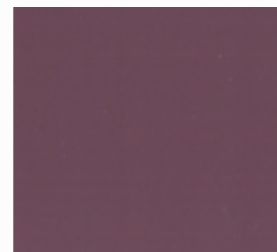

PLUM

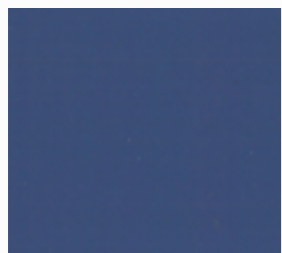

DEEP OCEAN



PETROLE



GREY

## COLOR CHART



### CLASSIFICATION

M1


Fine grain, mat and gentle to the touch.  
Rot-proof cotton base.

Composition : Mass-processed PVC 86%, Cotton 14%.

Thickness : 1,04 mm +/- 10%



Effective against  
bacteria(MRSA)

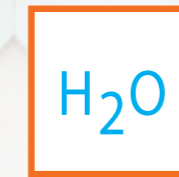

Solvent-free



Phtalate-free



Martindale abrasion  
1000000 cycles



Waterproof

A medical professional must always rely on his or her own professional clinical judgment when deciding whether to use a product(s) when treating a particular patient.

A healthcare professional anyway, refers to the package insert, product label and/or instructions for use before using any BMB MEDICAL product.

The products depicted are CE marked in accordance with EU Regulations and Directives.

BMB MEDICAL – 100 Avenue du Formans – CS 11001 – 01604 TREVOUX CEDEX – FRANCE

Tel : +33 474 08 71 71 – Fax : +33 474 08 71 72 – email : info@bmb-medical.com

Sarl au capital de 600 000 € - RCS Bourg en Bresse 89B652

SIRET : 950 413 229 00036 – Code NAF : 3250A – Code TVA : FR68 950 413 229

[www.bmb-medical.com](http://www.bmb-medical.com)



Tél : +33 (0)4 74 08 71 71

Fax : +33 (0)4 74 08 71 72

[info@bmb-medical.com](mailto:info@bmb-medical.com)

100 Av. du Formans - CS 11001

01604 TREVOUX Cedex

FRANCE

[www.bmb-medical.com](http://www.bmb-medical.com)